







# Epicutan-Testpflaster

Almirall Hermal GmbH

## TROLAB® Patch Test

Epicutan-Testpflaster  
Adhesive Strips for Patch Testing  
Bande adhésive pour tests épicutanés  
Cerotti adesivi per test epicutanei

Hersteller / Manufacturer / Fabricant / Produttore:  
Almirall Hermal GmbH  
Scholtzstraße 3  
D-21465 Reinbek  
Germany / Allemagne / Germania  
Phone: +49 (0)40 72704-0  
Fax: +49 (0)40 72704-317  
www.almirall.de

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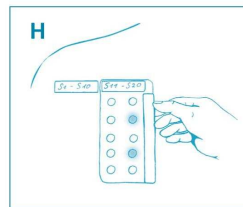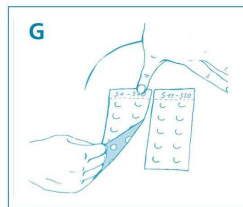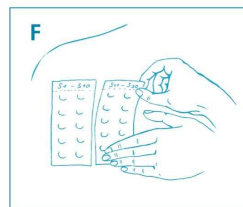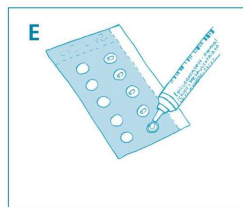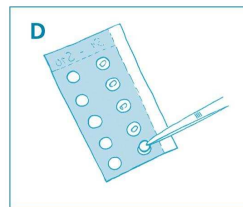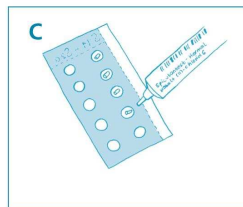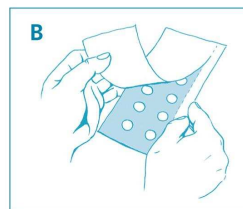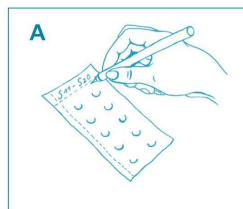#### Stand der Information

März 2011

## TROLAB® Patch Test Adhesive Strips for Patch Testing

1. Composition
2. Instructions for Use
3. Adverse Effects
4. Storage Instructions
5. Further Information

### 1. Composition

TROLAB® Patch Test is an aluminium-chambered test patch which is used as testing for contact allergies. This enables optimum occlusion during testing in accordance with the recommendations of the international professional associations.

Each test contains 10 aluminium chambers (internal diameter 8 mm, area 50 mm², volume approx. 20 µl) attached to a non-woven material consisting of rayon fibers. The adhesive used is an acrylate adhesive developed specifically for medical purposes.

### 2. Instructions for Use

#### Preparation

- Take a test patch and inscribe it on the top edge with the allergens to be tested (**Fig. A**).
- Remove the protective film from the patch – with the exception of the narrow strip of film on the right-hand side (**Fig. B**).

#### Filling the Test Chambers with Test Allergens in Petrolatum:

Place an approximately 5 mm-long quantity of test allergen ointment in Petrolatum directly – i.e. **without** filter paper – in the aluminium chambers (**Fig. C**). It is essential to avoid under-filling or overfilling the test chambers.

#### Filling the Test Chambers with Liquid Test Allergens:

- First you should attach a filter paper disc with a small amount of Petrolatum in the test chamber (**Fig. D**).
- Drip the liquid test substance (**Fig. E**) on to the filter paper immediately before testing (in order to avoid evaporation). Take up any excess liquid with absorbent material.

#### Adhering the Test Patch

- As a rule the filled patch is adhered to the back at a distance of 2-4 cm from the central line.
- Attach the lower edge of the test patch first. Spread it out towards the top with your hand so that it is fixed firmly in place and does not contain any air bubbles (**Fig. F**).
- Then press the individual test chambers firmly against the skin 2 to 3 times.

#### Removing the Test Patch

Pull the patch off below the perforation (**Fig. G**). The inscribed patch section on the upper edge of the patch can remain on the skin for the duration of the readings.

#### Reading

- If a **reading plate** is used, place the marking line of the plate on the lower edge of the inscribed patch section on the skin (**Fig. H**).
- If you read the test results **without a plate**, mark the areas with positive test reactions with a skin marker. This marking has to be renewed every time the results are read off.

### 3. Adverse Effects

- Allergic reactions to aluminium are very rare, although a few cases of contact allergies to aluminium have been described after hyposensitization injections with aluminium-adsorbed allergen preparations.
- When the patch is removed, the mechanical irritation may result in erythema with itching.

### 4. Storage Instructions

- Store at room temperature (10 - 27°C).
- Store in a dry place.

### 5. Further Information

#### Special Instructions

- Aqueous mercury preparations are incompatible with aluminium.
- Aluminium can promote the polymerization of acrylate monomers.

For the testing of these substances you should use test patches which do not contain aluminium chambers, e.g. Curatest®.

#### Pack Sizes

TROLAB® Patch Test is available in the following packs:

100 x 10 test chambers, with filter paper discs

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